
REF

ZORRO

Timago
International Group

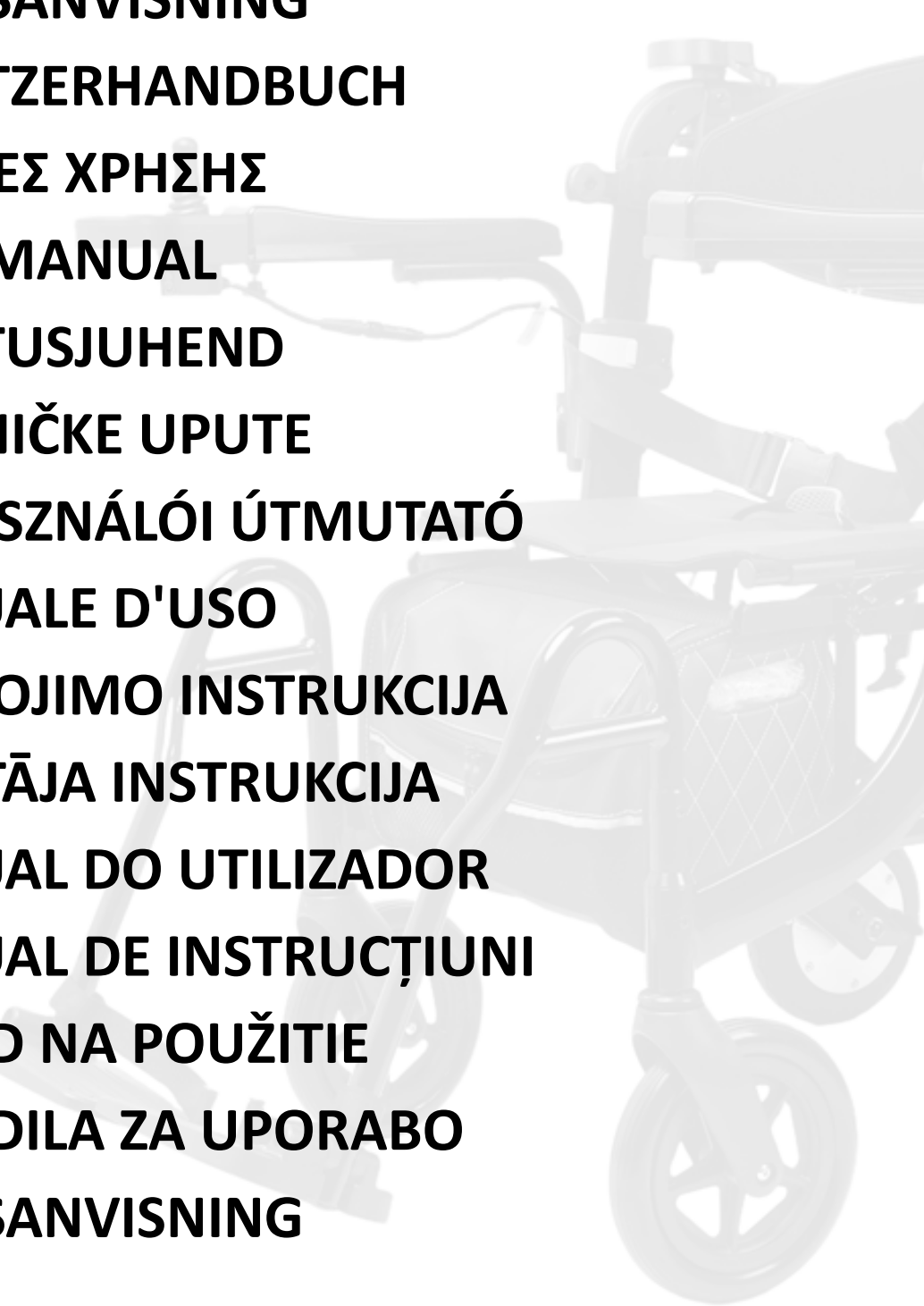PL	INSTRUKCJA OBSŁUGI
CS	NÁVOD K OBSLUZE
DA	BRUGSANVISNING
DE	BENUTZERHANDBUCH
EL	ΟΔΗΓΙΕΣ ΧΡΗΣΗΣ
EN	USER MANUAL
ET	KASUTUSJUHEND
HR	KORSNIČKE UPUTE
HU	FELHASZNÁLÓI ÚTMUTATÓ
IT	MANUALE D'USO
LT	NAUDOJIMO INSTRUKCIJA
LV	LIETOTĀJA INSTRUKCIJA
PT	MANUAL DO UTILIZADOR
RO	MANUAL DE INSTRUCTIUNI
SK	NÁVOD NA POUŽITIE
SL	NAVODILA ZA UPORABO
SV	BRUKSANVISNING

PL INSTRUKCJA OBSŁUGI	4
CS NÁVOD K OBSLUZE	18
DA BRUGSANVISNING	32
DE BENUTZERHANDBUCH	46
EN USER MANUAL	75
ET KASUTUSJUHEND	89
HR KORISNIČKE UPUTE	103
HU FELHASZNÁLÓI ÚTMUTATÓ	117
IT MANUALE D'USO	131
LT NAUDOJIMO INSTRUKCIJA	145
LV LIETOTĀJA INSTRUKCIJA	159
PT MANUAL DO UTILIZADOR	173
RO MANUAL DE INSTRUCTIUNI	187
SK NÁVOD NA POUŽITIE	201
SL NAVODILA ZA UPORABO	215
SV BRUKSANVISNING	229

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1. PRODUCT DESCRIPTION AND APPLICATION

Electric wheelchair is designed for people who have problems with mobility, helping disabled patients with their everyday tasks. It can replace a standard wheelchair.

Electric wheelchair is suitable for patients with:

- disability, who have problems with unassisted movement,
- dysfunctions and reduced movement,
- heart and cardiovascular system diseases,
- weakening body due to advanced age.

2. SAFETY INFORMATION

- Read the manual before first use.
- Do not use the wheelchair until fully assembled.
- Use of the wheelchair by people who are unable to operate it independently is not recommended.
- Maximum load capacity must not exceed 120 kg.
- The wheelchair can only be used by one person at a time.
- It is recommended to keep your feet on the footrest and your forearms on the armrests when driving.
- Always check the battery level on the control panel before use.
- This electric wheelchair should only be used on a flat surface or at low elevations. It should not be ridden on a surface inclined at an angle of more than 6 degrees.
- Do not drive the wheelchair at maximum speed when going uphill or downhill.
- Do not use the wheelchair on soft and slippery surfaces such as mud, sand etc.
- Do not ride the wheelchair sideways on highly elevated areas.
- Do not put any additional load such as a trailer etc. on the wheelchair.
- Do not repair or modify the construction of the wheelchair by yourself. Use only original spare parts provided by the manufacturer.
- Do not stand up or sit down on the wheelchair when the control module is switched on.
- When the freewheel (manual) mode is on, do not get on or off the wheelchair without assistance.
- Do not lean forward to grab items laying in front, on the side or behind the wheelchair. Do not stand on the footrest. It might result in the wheelchair falling over.
- Avoid rapid changes of the direction while driving on elevations.
- Always slow down while driving downhill or making turns.
- Before use check if the wheels are correctly mounted.
- Do not make sudden moves while using the joystick.
- When storing, protect the wheelchair from rain and adverse weather conditions.

- Avoid areas with horizontal gaps.
- If the wheelchair is not going to be used for a long period of time – remove the battery.
- Do not use the wheelchair on the road. It might result in a car accident.
- Do not store the wheelchair in places with high temperatures and humidity.

MEDICAL DEVICE INCIDENTS

Any serious incident involving a medical device must be reported to the manufacturer and to the competent authority of the European Union member state where the user or patient resides.

Hereby, we declare that the importer is not liable for any accidents, injuries, or damages resulting from prohibited actions stated in manual instruction. Damages caused by prohibited actions are not covered by warranty.

3. GENERAL INFORMATION

Do not operate wheelchair without first reading manual instruction and understanding it. Manual instruction provides user with necessary information regarding safety regulations, maintenance and operation of the vehicle.

Do NOT modify, add, remove or dismantle any of the wheelchair parts or disable functions by yourself. Contact authorized importer service for repairs.

PRE-RIDE SAFETY CHECK

Every time before using the wheelchair, perform inspection of the wearing parts such as breaks, wheels etc.

WEIGHT LIMITATIONS

Do NOT exceed weight capacity (120 kg) of wheelchair as it may result in damages that are no covered by the warranty.

Do NOT carry other passengers on the wheelchair, it may result in instability of a wheelchair or/and accident.

DRIVING ON INCLINES

Do NOT drive wheelchair across an incline or diagonally up or down an incline.

Do NOT drive up or down a potentially hazardous incline (e.g., areas covered with snow, ice, cut grass, or wet leaves).

DRIVING ON DIFFERENT SURFACES

Wheelchair construction allows driving on both outside and inside type of surfaces (concrete, asphalt, floor etc.). However it is important to remember about few rules to avoid unnecessary risk of damaging the wheelchair.

- Reduce speed while driving on uneven surfaces,
- Avoid driving on grass, mud or loosely packed gravel or sand,

- Do NOT drive on slippery surfaces such as ice as it may result in accident,
- Do NOT drive on surfaces that are unstable,
- Avoid holes and ditches,
- Do not drive on public roads as it may result in vehicle accident. Drive only on sidewalks or special lanes for bikers.

STAIRS AND ESCALATORS

Wheelchair is not designed to travel up or down the stairs or escalator. It can be transported via elevator.

PRECAUTIONS REGARDING WEATHER CONDITIONS

Avoid exposing the wheelchair to inclement weather conditions such as extreme heat or cold, rain, or snow. If you are caught in such conditions while riding, seek shelter and allow the wheelchair to dry completely before continuing or before storing or charging it.

Surface temperature can increase when exposed to external source of heat (e.g. sunlight).

BATTERY PROTECTION

Always protect the battery from freezing. **Do NOT** charge battery that has been frozen. If battery is damaged or cracked, immediately contact the importer service and do not use the battery.

Do NOT change any circuit connections. Any self-repairs are prohibited and might result damages that are not covered by the warranty.

SYMBOL	MEANING
	Safety Classification type B
	Warning! Read the manual instruction
	Prohibited! Performing a prohibited action can cause personal injury and/or equipment damage. Damages caused by prohibited actions are not covered by warranty.
	Mandatory! These actions must be performed before usage. If user fails to perform the mandatory actions, it may result in personal injury or/and equipment damage
IPX4	Protection against water spray from any direction
	Note / Read the manual instruction before operating

	Watch carefully for pedestrians
	Do NOT remove any parts or warning labels of the wheelchair
	Do NOT put fingers or limbs under covers or between movable parts. It may result in injuries
	Do not overload wheelchair. The maximum capacity of the scooter cannot exceed 120 kg
	Do NOT use wheelchair on stairs
	Do not use wheelchair in extreme weather conditions such as heavy rain, extreme humidity, icing etc.
	Warning! Wheelchair surface temperatures may be high
	Use only authorized importer service for repair in case of accident or/and damage to the parts
	Battery contains caustic substances. Dispose according to the local environment regulations
	Before driving, ensure that backrest is in upright position
	Use only lithium maintenance-free sealed batteries from manufacturer
	Keep tools and other metal objects away from the battery terminals. Contact with metal objects may cause an electrical shock

	Disposal and recycling Utilize according to local government laws.
	
	The battery charger is for indoor use only
	Flammable materials Do not expose to the open flame
	MR Unsafe
	Consult professionals before using wheelchair in case of mental or physical impairments
	Maintain proper position while driving
	Maintain your focus while driving and avoid distractions
	Fragile
	Protect from humidity and moisture

4. PRODUCT INFORMATION

A. ELEMENTS OF THE WHEELCHAIR

B. PRODUCT CHARACTERISTICS

- Aluminium, foldable frame reinforced with double crossbar
- Solid rear 8" and front, swivel 8" wheels
- Control module with LCD display
- Attendant module with power assist
- Safe and reliable electromagnetic brakes
- Attendant brakes with parking function
- Fully closed, maintenance-free battery
- Security belt and foldable sideway footrests
- Flip-up armrests
- Bag with hidden compartments for battery and engine
- Reflective lights on handles
- Adjustable backrest with 180° turn option
- Height adjustment for handles and backrest

C. TECHNICAL PARAMETERS

Wheelchair cannot be used in areas with flammable gas mixed with air, oxygen or nitrous oxide.


Note! Some parameters might vary depending on user weight, road conditions and battery usage.

Product code	ZORRO
Type	for indoor and outdoor use
Overall height	81-90 cm
Overall width	65 cm
Overall length	101 cm
Seat height from the ground	52 cm
Seat dimensions	45 x 35 cm
Folded dimensions	86 x 37,5 x 83 cm
Maximum capacity	120 kg
Battery weight	1,50 kg
Battery	lithium-ion, 24 [V] x 10 [Ah]
Battery charging time	4-6 h
Wheelchair weight	17,10 kg
Wheels	front: solid, swivel 8"

	rear: solid 8"
Maximum speed	≤ 6 [km/h]
Braking distance on horizontal surface	$\leq 1,5$ [m]
Maximum range	20 [km]
Max. angle of elevation	$\leq 6^\circ$
Min. turning radius	1,2 [m]
Safety classification	type B
Ingress protection level for engine	IPX4
Engine power level	2 x 250 [W]

D. ASSEMBLY AND ADJUSTMENT

- **Unfolding:** Remove the wheelchair from the box. Grab the handle located at the seat and slightly pull to unlock the folding lock (picture 1). Press seat on both sides till you hear click sound. Flip backrest to desired position (front/back). Install joystick frame into the opening under the armrest. Secure it with screw using hex key. Pick up desired length and tighten the securing screw afterwards (picture 2). Grab footrest frame and pull up, moving them inwards till footrest will be placed in front of the seat.
- **Folding:** To fold the wheelchair, grab by the handle on the seat and pull up. Press both sides of the wheelchair together till you hear click sound of the folding lock under the seat.
- **Control module adjustment:** Adjustments should be made while seated. Loosen the mounting screw located under the armrest, adjust the position or length of the module, and then tighten the screw afterwards.
- **Armrest tilting:** Press the clamp located on the outside of the wheelchair to lift the armrest. To lower the armrest, move it downwards till you hear click sound.
- **Backrest adjustment:** Backrest can be length adjusted and flipped by 180°. Unfasten the Velcro behind the backrest to adjust it length.

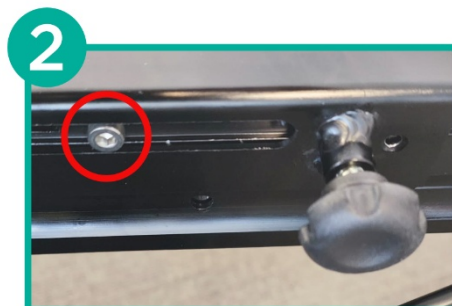5. OPERATION

A. MAIN CONTROL MODULE/ASSIST MODULE

- 1. Power button
- 2. Speed level (down)
- 3. Horn
- 4. Speed level (up)
- 5. Joystick
- 6. LCD Display (shows speed level and battery percentage)

- A. Resistance mode button
- B. Reverse mode
- C. Power button
- D. Power assist button
- E. LCD Display

B. ASSIST CONTROL MODULE FUNCTIONS

Wheelchair assist module applies automatically assist during ride, thanks to motion sensor. For more consistent power assist, press power assist button [D].

Display	Function	How to use function
-	Power button (turn ON/OFF button)	Hold and press power button to turn ON/OFF the assist module. While module is turned ON, pressing power button quickly again will change display information (ride time, total drive distance or current drive distance etc.).
	Brakes lock (parking mode)	When assist control module doesn't detect any movement of the hand within 15 cm range for 5 seconds, wheelchair will automatically change into parking mode (locking brakes). To leave parking mode, press any button on module.
  	Downhill mode 3 levels [M]	Downhill mode allows user to move wheelchair down the hill with higher resistance, helping prevent any accidents such as fall of the user from the wheelchair. To enter downhill mode press resistance mode button [A] and pick one of the three downhill speed options: 01, 02 or 03. 01 – slowest 03 – fastest The smaller the hill, the lower level of the mode should be used. Downhill mode can be used in reverse mode as well. Press reverse mode button [B] to enter reverse mode. For consistent apply of the power assist during ride, press and keep holding power assist button [D].
 	Uphill mode 3 levels [M]	Uphill mode allows user to move wheelchair up the hill with power assistance, helping assisting person push the wheelchair. To enter uphill mode press resistance mode button [A] and pick one of the three downhill speed options: C1, C2 or C3. C1 – slowest C3 – fastest The smaller the hill, the lower level of the mode should be used.

		Uphill mode can be used in reverse mode as well. Press reverse mode button [B] to enter reverse mode. For consistent apply of the power assist during ride, press and keep holding power assist button [D].
NON	Power assist [D]	To enter consistent assist mode, press and hold (not letting go) power assist button [D]. The display will show “NON” indicating that mode is active.
START	Motion sensor auto-assist	When assist module is turned ON, motion sensor will automatically sense assisting person movement within 15 cm. While grabbing handle in front of module, wheelchair will automatically enter auto-assist mode (“START” will show on display)
-	Reverse mode	Press reverse mode button [B] to use auto-power assist and modes while moving backwards.

C. CHARGER

- Set includes a charger dedicated to the battery.
- Battery charging time is around 4-6 hours.
- Plug the charger into the charging port located in the control module.

NOTE: Before plugging the charger in, turn off the control panel. The battery can also be charged directly. To do this, open the bag and remove the battery from the hidden compartment, then connect the charger to the charging port.

- Red light indicates that the battery is being charged, while green light informs that it is fully charged. After the charging is finished, the charger turns off to avoid overcharging.
- Connect the charger to power supply.
- Check the connection.
- To avoid the risk of fire, do not modify or disconnect the charging circuit while charging the battery.

D. DIRECTIONS FOR USE

- Tilt the footrest and sit comfortably on the seat. Lower the footrest back to horizontal position.

Note! Do not use the footrest as a support when sitting down or standing up from the wheelchair. This may cause the wheelchair to tilt and fall.

- Sit on the wheelchair and turn the control panel on. The LCD panel will light up. Joystick should be in neutral position.
- Take caution when using for the first time. Control the speed and direction carefully. The electromagnetic brake will get released only after pushing the joystick (electric

mode). To speed up, push the joystick further away from its neutral position. To slow down – do the opposite.

- To stop the wheelchair, release the joystick. Pressing the power button will also stop the wheelchair.
- Using dedicated buttons on the control panel, adjust the wheelchair's speed according to user's physical state and road conditions.
- Wheelchair should only be used on flat surfaces. Driving on muddy or uneven surfaces can cause control system damage.

6. MAINTENANCE INFORMATION

1. Wheelchair cleaning and maintenance

- It is recommended to regularly inspect the construction and its elements (wheels, frames, cables, control module, brakes etc.).
- Any used or damaged wheelchair parts need to be replaced (in such case, please contact the seller).

2. Battery maintenance:

- 1) Regularly check the battery level displayed. If the level is very low, charge it immediately to avoid a voltage drop (which may affect battery life).
- 2) If wheelchair won't be used for longer period of time, charge battery at least once a month.
- 3) Battery has a certain lifespan. Exchange your battery if it wears out.
- 4) Battery should be kept clean and away from kids. Avoid causing any shock or damage to the battery.
- 5) Charge it regularly to expand its lifespan.
- 6) Used batteries should be utilized according to local environmental regulations.

3. Storage and transport conditions

A. TRANSPORTATION

- Take caution during transportation. Do not throw, turn over or put pressure on the product and its packaging.
- For wheelchair folding, check „Assembly and adjustment”.

B. STORAGE

Wheelchair should be kept in dry and well-ventilated space. Avoid: high temperature and its rapid changes, acids, alkalis and corrosive chemicals.

C. STORING CONDITIONS

Temperature: -10[°C] ~+45[°C]

Humidity: ≤93%

Air pressure: 86 ~ 106 KPa

4. Troubleshooting

- In the event of a malfunction, the control module will indicate warning appearing on the screen. The number displayed on the screen will help to identify the problem.

Issue	Cause	Possible solution
E0002 E0004	E0002: Left motor malfunction E0004: Right motor malfunction	1. Check whether left or right motor main cables or any ports are damaged. If yes, contact seller service for their replacement. 2. Check if other cables are correctly connected. If all cables are correctly connected and are not damaged, contact seller service for further cause investigation.
E0008 E0010 E0020 E0040 E0800 E-021	E0008: Left motor detection error E0010: Right motor detection error E0020: Battery power level issues E0040: Malfunction (voltage) E0800: Control module error E-021: Rear control module error	Driver error. Contact seller service for repair.
EE	Emergency electric brake	Emergency braking happens when wheelchair speed exceeds the norm. Wheelchair will automatically lock motors for few seconds to stop wheelchair. Wheelchair will unlock motors automatically after sound alarm ends.
E-023	Issues with motion sensor	Contact the seller service for replacement of the rear control module.
E-024	Wheelchair calibration	After rear controller replacement, wheelchair might not be calibrated correctly. To calibrate wheelchair, press and hold power button on rear control module.
E-025	Motion sensor detection	Contact seller service to replace rear controller.
E-026	Motion sensor detection calibration	Contact seller service.
E-027	Display LCD issue	Contact seller service.

5. Warranty information

- All products distributed by our company are covered by a warranty, the terms of which are described in the warranty card available on our website. Please contact the retailer from whom you purchased the product. Please, keep in mind that for warranty purposes, the proof of purchase (receipt or invoice) should be retained. As a user-friendly company, we provide only tested products – in terms of material, quality and functionality. If you have any questions regarding service procedures, please contact us.
- **Warranty does not cover:**
 1. **wearable parts:** upholstery, tires, armrests, battery
 2. **damages resulting from:** self-repairs, changes in product construction, dropping, transportation, inappropriate usage, failure to follow instructions, natural disasters (fire, earthquake, flood).

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ZHEJIANG ZUBU MEDICAL TECHNOLOGY CO., LTD
SUNGO Cert GmbH